


Biotronix Healthcare®

BIOP80 Disposable Semi Automatic Biopsy System.

BIOP80 with detachable cannula, built-in Safety and detachable spacer.

BIOP80 represents the latest generation of semi-automatic precision cut biopsy needles

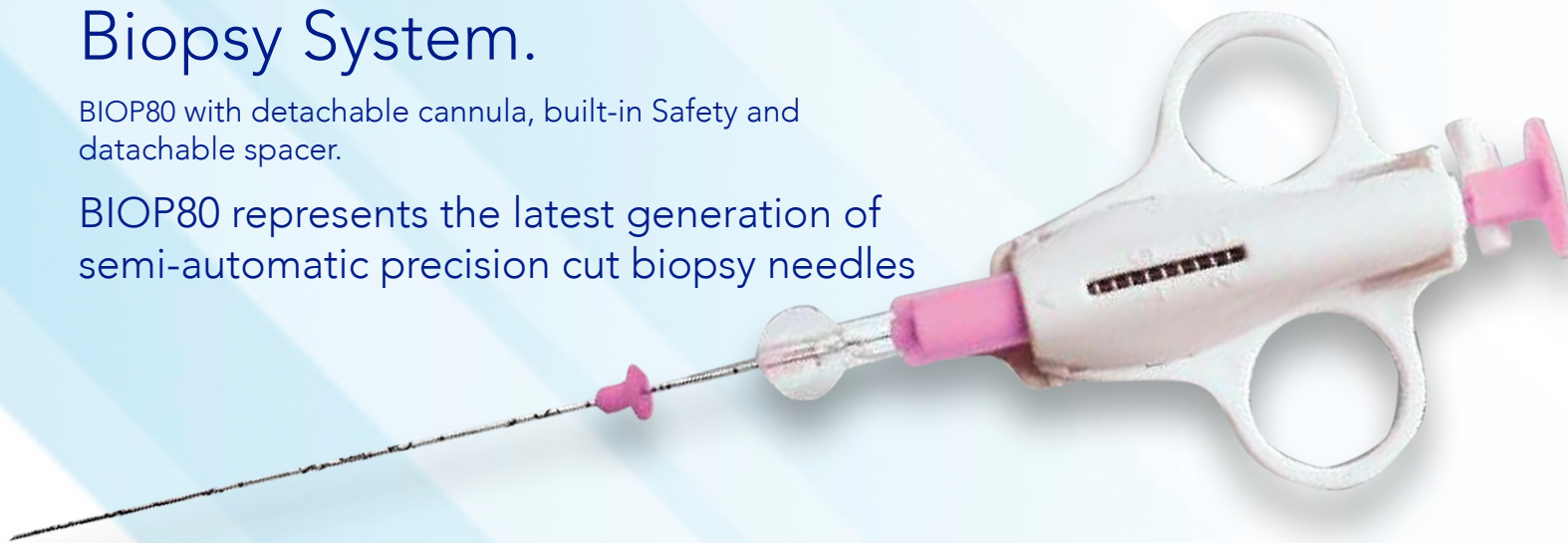

- **Unique Safety Lock.** BIOP80 is the first semi-automatic biopsy needles with built-in safety lock in the market which can lock the cannula to prevent accidental firing. This mechanism provides extra safety to patient during biopsy procedure, and operator when removing the samples from the needles.
- **Detachable Cannula:** BIOP80 is with detachable cannula. When used in CT guided procedure the main handle part of the needle won't obstruct in CT dome. The detached cannula and the stylet with pyramidal tip (packed within the same package) can be assembled as a coaxial needle which allows to perform multiple biopsies at the same puncture site.
- **Detachable Spacer:** The detachable spacer can prevent the needle from being "over penetrated". The spacer can be easily detached to further reduce the weight of the needle.
- **More Cutting Power:** The enhanced spring force and optimized internal design gives the needles 15%-20% more cutting power vs. other semi-automatic biopsy needles in the market. A powerful cutting force is even more important with the 10mm depth.
- **Two Penetration Depth:** 10mm and 20mm lengths can be selected adapting to the lesion size.

PRODUCT	GAUGE	LENGTH 10CM	LENGTH 16CM	LENGTH 20CM	PACKAGING
BIOP80 SEMI AUTOMATIC BIOPSY SYSTEM	14	BIO81028	BIO81032	N/A	10/BX
	16	BIO81029	BIO81033	BIO81037	10/BX
	18	BIO81030	BIO81034	BIO81038	10/BX
	20	N/A	BIO81035	N/A	10/BX

Manufactured For:
Biotronix Healthcare Industries Inc.
Miami, Florida USA